


Improve space-efficiency and profitability with the XYZ Distributor & Accuramatic Peristaltic Dispenser, for

Pharmaceutical Companies and Laboratories





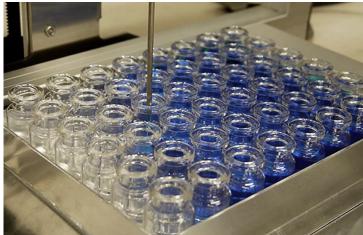

Industry-leading accuracy, to increase ROI

- Outstanding filling accuracy of \pm 0.5% improves your profitability by 50% (compared to the industry average of \pm 1%).
- Manual filling gives an average accuracy of \pm 10% overfilling just ten vials could mean giving away one full measure of product for free! We can prevent this for you.

Exceptionally hygienic, for international pharma-quality compliance

- Your product only comes into contact with the inside of one uninterrupted length of silicone tubing platinum cured, which will withstand repeated sterilisation.
- Some customers place the XYZ Distributor into a glove box or laminar flow box, to ensure a higher level of hygienic filling.
- Bacteriological inline filtration can be performed at the outlet side of the pump.

Space-efficient, for improved productivity

- Small footprint 790 x 625 x 390mm, to conveniently fit on a small workbench.
- Manually filling the same amount of containers could take 2 or 3 colleagues, reducing the cost-efficiency of the same production space by up to 66%.

Watch a video of the Accuramatic & XYZ >>>

Download the Accuramatic & XYZ Brochure >>>




Time-efficient, with fast changeover to fill many different products

- Fills an average of 12 15 containers per minute (for a 100ml fill).
- Can be operated by hand, foot switch, or be left to cycle automatically, leaving your operator free to assist production elsewhere.
- The disposable tubing can easily be removed and replaced with fresh sterile tubing, and the machine is instantly ready to fill your next product.
- Different vial sizes can be pre-programed and to re-called any time in the future. This significantly reduces the changeover time and also gives greater flexibility using the machine for different vial sizes.

Proven concept, over 400 units supplied worldwide

A high quality solution, appreciated by the most demanding customers, including Sigma-Aldrich, Allergan, Elemis, Schott, Teva Pharmaceuticals, NHS Hospitals and more.

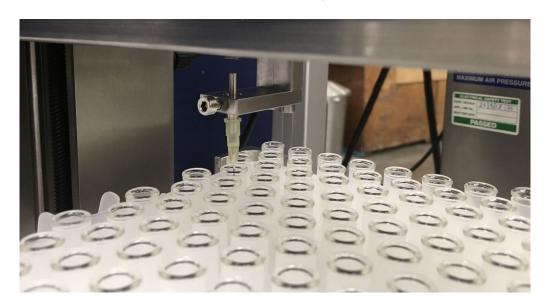

Our technical team will be delighted to suggest the best specification for you. Please share with us the below information, so that we can prepare a quotation:

- 1. Product Details
- 2. Product Viscosity
- 3. The container volumes required to be filled
- 4. The filling speeds expected

Contact our Sales Team for a quotation >>>

Important Notes:

Suitable for bottles and vials with a minimum neck diameter of 6mm. The products must be compatible with silicon tubing platinum cured. Not suitable for filling injectables.